HYPERFLUX

Operating Instructions

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SAFETY INSTRUCTIONS

Warnings

Read the entire instruction manual before using the Hyperflux Air Compression System.

No modification of this equipment is allowed.

If you experience severe pain, any unusual symptoms, or want to remove the attachments in an emergency during use:

- Stop the control unit by pressing the power button.
- Disconnect the hose from either the control unit or the attachments.
- Remove the attachments from your limbs.
- Promptly consult your licensed healthcare practitioner, as required.

Caution

- Do not attempt to take apart the system. The system has no user-serviceable parts. When service or repairs are required, please contact customer service at Hyperice.com
- Only use the charger provided with the system. Using a different charger may cause the system to operate incorrectly.
- To avoid risk of electric shock, do not use the system near water, such as near a bathtub, kitchen sink, laundry tub, or swimming pool.
- To avoid damage and risk of electric shock, never spill liquid of any kind on the system.
- Do not place the system, charger, or any accessories where they could be damaged, present a fall hazard, or become an obstruction to others.
- Keep the open ports of the control unit, hose interconnect, and power inlet free of debris.
- If the charger is damaged, the control unit is dropped or damaged, liquid is spilled on the system, or if the system does not operate normally when the operating instructions are followed, turn the system off by pushing the control unit's power button and then unplugging the system from the wall outlet. Contact customer service for assistance.
- Do not puncture or otherwise damage the attachments) since this may cause the system to operate incorrectly.
- To avoid risk of strangulation, do not leave a baby or child unattended with the power supply or hose.
- Choking hazard, small parts. Keep away from small children.
- Do not leave the system, charger, or any accessories where they could be damaged by children, pets, pests, or liquids. If you suspect your control unit is damaged, contact customer service for assistance.

- Do not allow lint or dust to accumulate on the control unit or the hose interconnects. If lint or dust accumulates, wipe down the system with a dry cloth before use.
- Do not walk while wearing the attachments.
- The attachments are designed to be used by only one person at a time.
- Do not hold the unit by the hose.
- Consult your physician before using this product if you are experiencing inflammation, an infection, pain of unknown origin, bleeding (internal or external) at or near the site of application, or if you have a wound at or near the site of application.
- Consult your physician before using this product on sensitive skin.
- Do not use this product without consulting your physician if you have any of the following:
 - Acute pulmonary edema
 - Acute thrombophlebitis
 - Acute congestive cardiac failure
 - Acute infections
 - Deep vein thrombosis (DVT)
 - Episodes of pulmonary embolism
 - Wounds, lesions, or tumors at or near the site of application
 - Where increased venous and lymphatic return is undesirable
 - Are pregnant, have diabetes, neuropathy, have bone fractures or dislocations at or near the site of application
- Do not use the Hyperflux air output hose to direct pressurized air toward your eyes, nose, mouth, or ears. Doing so may lead to serious injury.
- Use by unconscious or incapacitated persons may be dangerous without supervision.
- Make sure the power inlet on the control unit is easily accessible at all times in order to disconnect power if required.

Save this manual for future reference.

ILLUSTRATIONS

Hyperflux Control Unit


1. Power button
2. Air outlet and power inlet
3. Time adjustment button (decrease)
4. Time adjustment button (increase)
5. Pressure adjustment button (decrease)
6. Pressure adjustment button (increase)
7. Display screen
8. Play/pause button

Hose


1. Junction box air outlets
2. Junction box
3. Blocking plug (on the underside of the junction box)
4. Connector

Leg Attachment

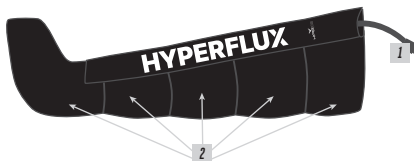1. Attachment connector
2. Zones

Charger

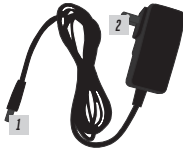1. Barrel connector

2. Wall outlet plug

OPERATING INSTRUCTIONS

Adjust the Intensity

Adjust the intensity level of the session by pressing the  and  buttons on the left and right of the level indicator. Intensity level 1 is the gentlest setting. The air compression becomes stronger as the intensity level is increased.

Intensity can be adjusted while the session is running. When level is adjusted during a session, the system will stop for 10 seconds, and the cycle will resume at the new pressure.

Adjust the Session Time

Adjust the session time by pressing the  and  buttons on the left and right of the time indicator. The session time can be set between 10 minutes and 2 hours and 55 minutes (in 5-minute increments). A typical session time is between 15 and 60 minutes long.

Setting the time longer than 2 hours 55 minutes will put the control unit into Continuous Mode. To turn off Continuous Mode, decrease the time. In Continuous Mode, the screen will display CONTINUOUS instead of the timer and will begin counting up from zero when the session starts.

Time can be adjusted while the session is running. Tap the time adjustment buttons to add or subtract time from the session in 5-minute increments.

Start the Session

To start the session, tap the play/pause button.

Stop or Pause the Session

To stop the session at any time, tap the play/pause button. This will pause the session. To restart your paused session, tap the play/pause button again.

Finish the Session

The session will continue massaging until time runs out and the display reads 0:00:00. The time indicator will blink, and the control unit will continue until the current cycle is finished.

When the session is completed, disconnect the hose from the control unit and the attachments, remove the attachments, and turn off the control unit by pressing the power button.

To disconnect connectors from the junction box or the control unit, push the button on the top of each connector while pulling away.

Turn Off the Control Unit

To turn off the system, press the power button and confirm that the green power LED is off.

Adjust the Rest Time

To enter rest time edit mode, set the time to 10 minutes then hold down the  time adjustment button for four seconds (until the red rest time indicator lights up and begins blinking).

When in rest time edit mode, the rest time will blink on the screen. Use the  and  time adjustment buttons to increase or decrease the rest time. The rest time can be set between 15 and 90 seconds in 15-second increments. The default rest time is 30 seconds.

To save and exit rest time edit mode, tap the play/pause button on the control unit.

Air Compression Patterns

The air compression pattern can be set to Pulse or Sequential. If the control unit is set to SEQUENTIAL mode, it will be displayed under the session time. To change the air compression pattern, set the control unit to Level 7 and then hold the  level button for five seconds.

Before the Pulse air compression pattern begins, you will experience a pre-inflate cycle, during which the connected attachments are molded to your exact body shape. Once the pre-inflate cycle is complete, the Pulse air compression pattern will begin by compressing your feet. Each zone of the attachments will first compress in a pulsing manner and then release as the air compression pattern works its way up your limb. When the top zone completes its air compression, there will be a brief rest period and then the cycle will begin again. This will repeat until the session time runs out. When a Pulse air compression session is resumed after a pause, the system will perform the pre-inflate cycle before continuing.

The Sequential air compression pattern begins by compressing the bottom zone. Each zone of the attachments will compress and hold pressure. This will repeat for each zone of the attachments as the air compression pattern works its way up your limb.

Adjust the Number of Zones

The number of attachment zones enabled can be changed between one and five zones. Zones can be disabled from the top to the bottom of the attachment, one zone at a time. Zone adjustments do not reset after powering off the control unit.

To enter zone edit mode, set the control unit to Level 1 and then hold the  level button for five seconds. The top-most enabled zone of the attachment will begin blinking.

When in zone edit mode, the top enabled zone on the screen's attachment graphic will blink. Use the  and  buttons on the control unit to increase or decrease the number of enabled zones. At least one zone must be enabled.

To save and exit zone edit mode, tap the play/pause button on the control unit. The top-most enabled zone will stop blinking.

CARING FOR THE SYSTEM

Cleaning the System

To clean the control unit:

- Wipe down the system with a damp, clean cloth.
- Dry thoroughly with a clean cloth.

Cleaning the attachments:

- Wipe down the attachments inside and out with a damp, clean cloth.
- Dry thoroughly with a clean cloth.
- Do not machine wash or dry.
- Do not dry clean.

Maintaining the System

The control unit, hose, charger, and attachments require no routine maintenance or service except for the care in this section.

Storing the System

Store control unit, hose, charger, and attachments in a clean, dry location.

Replacement Parts

Please contact customer service at Hyperice.com.

TECHNICAL INFORMATION

Do not attempt to take apart the system.

The system has no user-serviceable parts. There are no user-replaceable fuses.

Internal Battery Information

This control unit is equipped with a rechargeable lithium ion battery. The internal battery is designed to allow use of the Hyperflux anywhere—even when power outlets aren't available. The Hyperflux may need to be plugged in before first use. The battery will provide power for 2+ hours of continuous use. It takes approximately 6 hours to fully charge the battery.

The rechargeable lithium ion battery is intended to be changed only by authorized service personnel.

Product Specifications

- Hyperflux Dimensions: 4" (width), 4.5" (depth), 8" (height); [10.2 cm (width), 11.43 cm (depth), 20.32 cm (height)]
- Hyperflux Weight: 3.4 lbs [1.54 kg]
- Hyperflux electrical requirement: 15V  DC 1A

AC-DC Charger

Only use the AC-DC charger model number PSAW30I-13-BNL24 provided with the system. Using a different charger may cause the system to not operate correctly.

- Input: 100-240V  0.8-0.4 A 50/60 Hz per Hyperice model number PSAW30I-13-BNL24
- Output: 15V  DC minimum 1.6 A per Hyperice model number PSAW30I-13-BNL24

TROUBLESHOOTING

Problem	Possible Causes	Solutions
The system does not start	Power is not turned on Charger is not securely connected Faulty electrical wall outlet	Press the power button to turn the control unit on. Check that the charger is securely connected to the Hyperflux control unit and the electrical outlet. Check that the wall outlet works.
The attachments do not inflate	The session has not been started The hose is not securely connected The attachments have been damaged	Tap the start button to start the session. Check that the hose is securely connected to the Hyperflux control unit and that the attachments are securely connected to the junction box. Check that the attachments do not have an air leak.
The system stopped pumping	The hose is not securely connected The attachments have been damaged	Check that the hose is securely connected to the Hyperflux control unit and that the attachments are securely connected to the junction box. Check that the attachments do not have an air leak.
Air leak message: Err	Air leak	Check for leaks in the hose or attachment. Check that the connectors are firmly connected.
Low Battery	Battery needs to be charged	Plug in the control unit to charge the battery.

Contact customer service at [Hyperice.com](https://www.hyperice.com).

WARRANTY

IF WITHIN ONE YEAR FROM THE DATE OF PURCHASE, THIS PRODUCT FAILS DUE TO A DEFECT IN MATERIAL OR WORKMANSHIP, HYPERICE, INC. WILL REPAIR OR REPLACE THE PRODUCT, OR NECESSARY COMPONENTS, FREE OF CHARGE.

This Warranty Excludes:

- (a) damage caused by accident, abuse, mishandling, or transport;
- (b) units subjected to unauthorized repair;
- (c) units not used in accordance with Hyperice care instructions;
- (d) damage exceeding the cost of the product;
- (e) deterioration of the delivered product resulting from abnormal storage and/or safeguarding conditions on the client's premises, and
- (f) failure to provide the dated proof of purchase.

Some states and countries do not allow a limitation of damages, so the foregoing limitation may not apply to you.

This warranty guarantees specific legal rights, and but other rights may vary from country to country and from state to state. This warranty service is available through the operating policies and procedures of Hyperice Inc. and all Hyperice international distributors. If you are a client, please contact the appropriate Hyperice business entity, dealer or retailer, from whom you directly purchased the product for warranty and return authorization procedures. This warranty is initiated and executed by Hyperice Inc.

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Copyright 2020 © Hyperice, Inc. All rights reserved.
Hyperice, Inc.